



HYDRIM[®] G4

Washer-Disinfector






LEADERSHIP

in Infection Control

Instrument washing is the first step in the sterilisation process, which minimizes the physical handling of instruments and the many steps involved with manual cleaning such as pre-soaking, scrubbing and drying instruments. In addition to meeting the requirements of various guidelines and controls, automated washing and disinfection provides a safe, effective and efficient way to clean instruments, allowing staff to focus on the most important thing in the dental office, the patient.

Through a sustained legacy of product innovation and refinement, SciCan developed the HYDR/M[®] Washer-Disinfectors, which provide a suitable solution for nearly every office, whether as a stand-alone or counter-top unit.



Originally founded in 1957, SciCan develops and manufactures highly innovative infection control products that address the need for speed, flexibility, and efficiency in today's environments.

In 1989, SciCan developed and produced an internationally recognized infection control innovation, the STAT/M® Cassette Autoclave. Over the next quarter century, SciCan continued to be at the forefront of technology with the release of the HYDRIM® Instrument Washer-Disinfectors, BRAVO™ Chamber Autoclaves, and G4 Networking Technology.

SciCan is a global leader in infection control specializing in dentistry and ophthalmology in more than 100 countries around the world. The portfolio of products continues to grow to address market needs; innovative products such as the SALUS® Instrument Reprocessing container, high-quality turbines, handpieces and contra angles with the appropriate cleaning and maintenance unit, the STATMATIC™ Smart, and the OPTIM® Cleaner and Disinfectant complete the product range.

In 2007, SciCan opened the European headquarters, located in Germany. Our team consists of highly motivated and technically informed individuals whose network of authorities, infection control specialists and other manufacturers allow us to provide legally compliant and evidence-based information. We work closely with specialists in the dental and medical field, enabling us to provide practical guidance and help drive product improvements that will benefit our customers.

It is very important to me and my team, especially in the time of online sales, to offer you not only the best products but also a service that makes our customers feel well advised and cared for. Please don't hesitate to contact our office for any infection control related questions.

Good luck. I wish you only good days and sufficient time for the essential things.

Stefan Helsing (Managing Director)

HYDRIM G4

Instrument Washer-Disinfectors

HYDRIM G4 automated instrument washer-disinfectors are the first step in the sterilisation process.



Reliable

Both models are designed in accordance with EN ISO 15883. The washer- disinfectors have validated cycles to ensure repeatable performance. HYDRIM meets all EU Standards and RKI Requirements related to washer-disinfectors.



Safe

Automated washer-disinfectors reduce physical handling of instruments, decreasing the risk of exposure and injuries. HYDRIM G4 ensures that every set of instruments receive the same validated cleaning process by dosing the correct amount of detergent with each cycle.



Dry

Active air drying technology with hepa filter quickly dries the instruments and allows them to immediately be transferred to the next step of instrument reprocessing, sterilisation.



Capacity

Two unit sizes (table-top and under-the-counter configurations) with different options for baskets and adaptor ensure you have a unit that meets the needs of any practice. HYDRIM M2 G4 has a large chamber capacity, allowing offices to turn over 1500 instruments in a 10 hour day.* HYDRIM C61WD has a 6 full-size cassette capacity, allowing offices to turn over 120 instruments in just 60 minutes – start to dry.*



Compact Design

The HYDRIM C61WD, as a table top unit, is perfect for smaller offices and the HYDRIM M2 G4 is specially designed for medium and large offices with 5 or more operators. It can be installed under the counter or as a freestanding unit. The cleaning solution is stored in the unit, so there is no extra storage space required.



G4 Technology

Allows for reliable load release; connects and sends cycle data to your smart devices or computer, providing easily accessible cycle data, maintenance information and service logs. Should you wish, G4 technology can connect your HYDRIM to specialized technicians, which allows for remote support.

*Actual cycle times will vary with incoming water temperature and pressure.

HYDRIM M2 G4



HYDRIM C61WD G4

Cleaning

Cleaning is one of the most important steps in the instrument reprocessing procedure. In order to achieve effective disinfection and sterilisation, the instruments need to properly be cleaned first. Washer-Disinfectors have been proven to be the most safe and effective way to clean and disinfect instruments, which is supported by many regulating authorities and dentists worldwide.



The powerful cleaning and disinfection units are available in two sizes, HYDRIM M2 G4 and HYDRIM C61WD. Both units are made of robust stainless steel and offer, with a comprehensive range of accessories, solutions for almost every office, either as table-top or under-the-counter unit. No extra detergent space is required as the patented cleaning solution is stored inside the unit. HYDRIM washer-disinfectors are user-friendly, safe and effective, allowing for thorough cleaning as well as thermal disinfection at 93°C.

By default all HYDRIM units offer active air drying with a hepa filter, protecting instruments from recontamination and corrosion. As standard, the trolleys and the dosing pump for the liquid cleaner are also included in the scope of delivery.

Conformity

HYDRIM M2 G4 and C61 WD G4 are produced in conformity with the following directives and standards:

| | | | |
|---------------------------|--|--------------------------|---|
| 93/42/EEC | Medical device class IIa (European Directive for Medical Devices) | EN 1717 | Protection against Pollution of Potable Water |
| EN ISO 15883-1/-2 | Washer-Disinfectors | EN ISO 13485:2016 | Medical devices – Quality management systems – Requirements for regulatory purposes |
| EN 61010-1 /-2-040 | Safety Requirements for Electrical Equipment for Measurement, Control and Laboratory Use | EN ISO 14971:2012 | Application of risk management to medical devices |
| EN 61326-1 | Electrical equipment for measurement, control and laboratory use – EMC requirements | KIWA certified | |

G4 Technology

As special feature of HYDR/M units is the innovative G4 technology. A large touch screen provides in-depth cycle info and allows users to view current operation.

With an integrated Ethernet port, HYDR/M can connect to the office network and send cycle data to your smart devices or computer, providing easily accessible cycle data, maintenance information and service logs.

To help optimize office processes, the HYDR/M's can be configured to email cycle information or error messages directly to office staff or service technicians. Cycle data is automatically saved internally, on a USB and can also be emailed or printed with a network printer, a feature that has been specifically designed for meeting the growing legal demand of physical documentation.

The HYDR/M G4 units have a traceability features that requires a user code when releasing a load.



HYDRIM M2 Specifications



| | | |
|-------------------------------------|---|--|
| Size | Dimensions: (w x h x d) | 598 x 850 x 600 mm |
| | Depth with door open | 1200 mm |
| | Weight | 90 kg |
| Cycle Times* | Rinse and Hold | 18 mins |
| | Regular Wash + Dry** | 35 mins. + 30 mins. |
| | Heavy Duty Wash with Disinfection + Dry | 60 mins. + 30 mins. |
| Water Consumption Per Cycle† | Rinse and Hold | 16 litres |
| | Regular Wash | 32 litres |
| | Heavy Duty Wash & Disinfection | 40 litres |
| Electrical Connection‡ | Electrical Rating | 230-240V, 50 Hz or 60 Hz, 12A or 13A, single phase |
| Water and Drain Connections | Hot and Cold Water Tap Connection | 3/4" Garden hose, male |
| | Recommended Hot Water Temperature | 60°C |
| | Cold Water Temperature | < 30°C |
| | Inlet Water Pressure | 2 - 5 bar |
| | Drain | 2 cm |

* Actual cycle times will vary with incoming water temperature and pressure.

** This cycle DOES NOT provide disinfection and is not available for use. This cycle can be added by a technician.

† Water consumption includes drying

‡ Electrical requirements are model specific

HYDRIM M2 Rack Configurations

Configuration 1



Vertical Instrument Rack
01-110411S



Hygiene Basket
01-109967S




Basket with Hinged Lid
01-109966S

Configuration 2



Vertical Instrument Rack
01-110411S

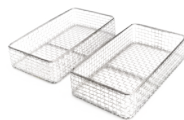

Mesh Basket
01-114818



Bar for Scaler Tip with Torque Tool
01-114841




LCS Adapter
01-114405S

Additional Accessories

Adapter for Handpieces
01-114845

Metal Plate Holder for Bars
01-114846

Configuration 3

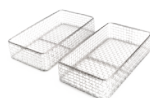

Mesh Basket
01-114818



STATIM 2000 Basket
01-107240



Bar for Scaler Tip with Torque Tool
01-114841



Multi-Use Bar Accessory
01-114842



LCS Adapter
01-114405S


Additional Accessories

Adapter for Handpieces
01-114845


Metal Plate Holder for Bars
01-114846

Adapter for Surgical Suction Tubes
01-114844

Configuration 4

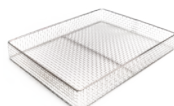

Mesh Basket
01-114817

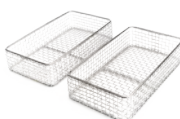

Mesh Basket
01-114818



LCS Adapter
01-114405S

Additional Accessories

Adapter for Handpieces
01-114845

Metal Plate Holder for Bars
01-114846

More loading samples available upon request.

HYDRIM C61WD Specifications



| | | |
|-------------------------------------|---|---|
| Size | Dimensions: (w x h x d) | 598 x 520 x 526 mm |
| | Depth with door open | 829 mm |
| | Weight | 44 kg |
| Cycle Times* | Rinse and Hold | 5 mins |
| | Regular Wash + Dry** | 26 mins. + 10 mins. |
| | Heavy Duty Wash with Disinfection + Dry | 42 mins. + 10 mins. |
| Water Consumption Per Cycle† | Rinse and Hold | 11 litres |
| | Regular Wash | 11 litres |
| | Heavy Duty Wash & Disinfection | 30 litres |
| Electrical Connection‡ | Electrical Rating | 208-240V or 200-230V, 50 Hz or 60 Hz, 12A, single phase |
| Water and Drain Connections | Hot and Cold Water Tap Connection | 3/4" Garden hose, male |
| | Recommended Hot Water Temperature | 60°C |
| | Cold Water Temperature | < 30°C |
| | Inlet Water Pressure | 1 - 10 bar |
| | Drain | 2 cm |

* Actual cycle times will vary with incoming water temperature and pressure.

** This cycle DOES NOT provide disinfection and is not available for use. This cycle can added by a technician.

† Water consumption includes drying

‡ Electrical requirements are model specific

HYDRIM C61WD Rack Configurations

Configuration 1



3 Cassette /
2 Basket Rack
01-113253



Handpiece
system LCS
01-113625



Vertical
Instrument Rack
01-113255



Bar for Scaler Tip
with Torque Tool
01-114841

Additional Accessories

Adapter for
Handpieces
01-114845

Metal Plate
Holder for Bars
01-114846

Configuration 2



Handpiece
system LCS
01-113625



Vertical
Instrument Rack
01-113255



Bar for Scaler Tip
with Torque Tool
01-114841



STATIM 2000
Basket
01-107240



Basket Rack
01-113257



Lower Trolley
Basket
01-114942

Additional Accessories

Adapter for
Handpieces
01-114845

Metal Plate
Holder for Bars
01-114846

Configuration 3



3 Cassette /
2 Basket Rack
01-113253



Handpiece
system LCS
01-113625



Bar for Scaler Tip
with Torque Tool
01-114841



Multi-Use Bar
Accessory
01-114842

Additional Accessories

Adapter for
Handpieces
01-114845

Metal Plate
Holder for Bars
01-114846

Adapter for Surgical
Suction Tubes
01-114844

Configuration 4




Full Basket
01-113545



Basket Rack
01-113257



STATIM 2000
Basket
01-107240

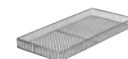

STATIM 5000
Basket
01-107241

More loading samples available upon request.

HYDROS™ - Water Treatment System

for use with washer-disinfectors

Bad water quality can lead to discoloration and calcium deposits in the machine. Connecting the washer-disinfector to a R/O water treatment system is therefore recommended to achieve optimal cleaning and support the longevity of your machine.

Demineralized water (R/O water) is also recommended in the RKI guideline from 2012 (Hygiene requirements for the reprocessing of medical devices) for the final rinse in washer-disinfectors.

HYDROS has a special multi-layered filter that removes several unwanted components of the water, step by step.



User Friendly Operation

Easy exchange of the filter cartridge



Space Saving Installation

Horizontal and vertical installations possible



Easy To Use

Suited for units with hot water connection
(up to 60°C)

Technical Specifications

| | |
|--|----------------|
| Capacity at a carbonate hardness of 10° dH * | 5000 l |
| Operating pressure | 6 bar |
| Water intake temperature | 4 – 60°C |
| Flow at 1 bar pressure loss | 850 l/h |
| Nominal flow | 300l/h |
| Pressure loss at nominal flow | 0.45 bar |
| Dimensions (width/depth/height) | 288/255/550 mm |
| Weight (dry/wet) | 18 / 24 kg |
| Connections (input/output) | G1" / G3/4" |

* The stated capacity is a guide value only, that may vary by +/- 20% depending on local water quality and machine type.

SciCan's Full Spectrum Solutions

Helping you get the most out of your sterilisation area and the equipment you put in it.





Manufactured by:

SciCan Ltd.

1440 Don Mills Road, Toronto,
ON, Canada M3B 3P9

Tel: 416-445-1600

Fax: 416-445-2727

Toll free: 1-800-667-7733

customerservice@scican.com



EU Representative & Distributed by:

SciCan GmbH

Wangener Strasse 78

88299 Leutkirch

Deutschland

Tel.: +49 (0)7561 98343 - 0

Fax: +49 (0)7561 98343 - 699

customerservice_eu@scican.com

Office for Swiss Customers:

SciCan Medtech AG

Alpenstrasse 14

6300 Zug

Switzerland

Tel: +41 (0) 41-727 7027

Fax: +41 (0) 41-727 7029