

INTRODUCING

THE REVOLUTIONARY
HANDPIECE ATTACHMENT
FOR MOCOM TETHYS H10+



THANKS TO THE NEW HMD, TETHYS H10+ IS NOW ALSO A HANDPIECE THERMAL DISINFECTOR.

The new Handpiece Attachment for Mocom Tethys H10+ allows handpieces and scaler tips to be thermally disinfected. This latest accessory makes it possible to wash, thermally disinfect and dry up to 6 instruments, turbines or micromotor handpieces, simultaneously. You can also process up to 3 scaler tips at once.

The Handpiece Attachment is compatible with all Mocom Tethys H10+ units, including those that have already been installed, a simple firmware update is all that is needed.



Attachment easily removed.

ADAPTERS FOR ANY NEED.

The Handpiece Attachment allows you to process:

- Micromotor handpieces
- Turbines
- Scaler tips

of instruments, in order to ensure proper cleaning of the instrument.

Just insert the specific adapter for the family

Tethys H10+With the new handpiece cycle and the dedicated support, a cycle can be

Handpiece Cycle on your Mocom

completed in just 34 minutes.

This cycle consists of:

 Washing internal & external (no ultrasound)

• Thermal disinfection

Drying











Phase 1 Washing with Detergent

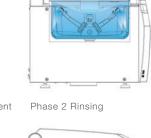



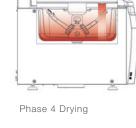

Phase 3 Thermal Disinfection

FOR MOCOM TETHYS H10+, CONTACT MOCOM AUSTRALIA FOR YOUR FIRMWARE UPDATE TODAY.

TO START USING THE HANDPIECE ATTACHMENT

0427 816 459 or sales@mocomaustralia.com.au with any questions.

mocomaustralia.com.au to book in your firmware upgrade or direct at

Please contact your local Mocom Australia Service Technician via



1300 65 88 22 | henryschein.com.au



At Henry Schein Dental, our mission is to improve the lives of all Australians by focusing on practice care, so dental professionals can focus on patient care.